



Reichert 7

Auto Tonometer

One-touch, non-contact tonometry.



Reichert 7, Reichert Technologies' seventh generation non-contact tonometer, is the most user friendly NCT ever. With a number of new features, including an easy to navigate touch screen interface and multiple-measurement mode, Reichert 7 will be a welcome addition to your practice. The totally automated operation results in error-free measurements and high patient throughput. Measurement speed and simplicity without sacrificing accuracy.



Multiple-measurement mode:

Reichert 7 features an all-new, one-touch, multiple-measurement mode that delivers 3 or 4 measurements with a single press. This new feature saves time, increases accuracy, and ensures the highest level of comfort for your patients.

User interface:

The touch screen interface features an easy-to-use, icon-based operating system enabling even inexperienced users to make measurements in seconds. The full color display guides the operator through the measurement process showing everything from the alignment of the eye, to the display of the readings. Set up menus and instrument options are also easily accessible and simple to navigate.

Soft puff and quiet operation:

Reichert 7 is designed to make the entire measurement process a comfortable and non-intimidating experience for patients. The instrument features the quietest operation and softest air puff of any NCT available.

Left/Right sliding headrest:

A sliding headrest enables easy patient positioning to ensure fast and comfortable measurements. The operator simply slides the headrest to the left or right position and the patient positions their head against the comfortable headrest pad. Reichert 7 automatically detects which eye is being measured.

Bright LED array for fast fixation:

The Reichert 7 nosepiece contains an array of bright LEDs. These provide an obvious visual cue to the patient for fast fixation and alignment, thereby reducing exam time.


Accuracy, simplicity, speed:

Operation of Reichert 7 is extremely quick. Patients simply lean against the forehead rest while the operator touches an icon on the screen to activate the fully automated alignment and IOP measurement process. The measurement is taken with speed and simplicity unlike any other NCT, without sacrificing accuracy.

¹Optional motorized chinrest is available (see Catalog Number 16049)




Headrest enables easy patient positioning

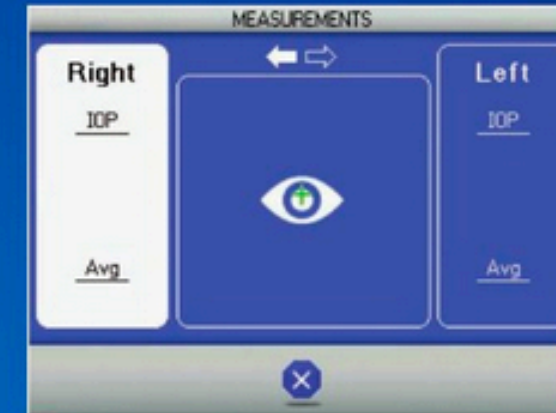

Auto-alignment positioning



No joystick, chinrest or elevation controls necessary¹

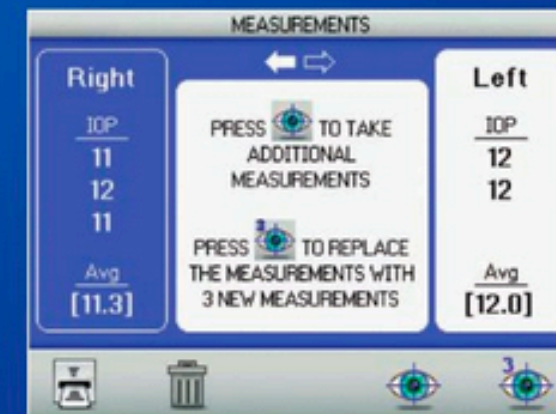

Clearly displayed measurement options and data



Features

- Alignment is completely automated, no need for joysticks
- No chinrest or elevation controls necessary – patients simply lean against the forehead rest
(Optional motorized chinrest is available)
- Easy-to-use, touch screen user interface
- One-touch, single or multiple-measurement mode
- Clearly displays all measurement data
- EMR connectivity via USB port
- Internal printer simplifies record keeping

Reichert 7 Specifications

Catalog Number	16050
Height	19.75 in (50.2 cm)
Width	10.5 in (26.7 cm)
Depth	14 in (35.6 cm)
Weight	23 lbs (10.4 kg)
Voltage	100/240 VAC
Frequency	50/60 Hz
Measurement Range	0 - 60 mmHg (7 - 60 mmHg ISO 8612 tonometer standard)



Optional Chinrest Specifications

Catalog Number	16049
Vertical Travel	65 mm, motorized

Proudly Made in the USA

Call your authorized Reichert Technologies Distributor today.

Important Notice

For your own protection, only purchase your Reichert Technologies instruments from and have them serviced by an Authorized Reichert Distributor.

16050-110-Rev-A ©2011 Reichert, Inc. Specifications subject to change without notice. Reichert is a registered trademark of Reichert, Inc. Reichert Technologies and "Advancing Eye Care. Preserving Sight." are trademarks of Reichert, Inc.



Reichert, Inc.
 3362 Walden Avenue, Depew, NY 14043
 Toll Free: 888-849-8955 • Tel: (716) 686-4500
 Fax: (716) 686-4545 • www.reichert.com

ISO 9001 Certified

Reichert products are designed and manufactured under quality processes meeting ISO 9001 requirements.